

Reference Material Certificate

741/01

Aluminium Base (Type of Standard) AlZn, Set 740

Certified Values

Element	Mass content [%]	Uncertainty [%]
Silicon (Si)	0.31	± 0.01
Iron (Fe)	0.40	± 0.02
Copper (Cu)	0.10	± 0.01
Manganese (Mn)	0.105	± 0.007
Magnesium (Mg)	0.035	± 0.003
Chromium (Cr)	0.023	± 0.002
Nickel (Ni)	0.011	± 0.001
Zinc (Zn)	0.805	± 0.025
Titanium (Ti)	0.027	± 0.002

The uncertainty reported is the result of standard deviation of all results multiplied with a factor of two and represents approximately the 95% confidence interval.

Manufacturing

This certified reference material for the analysis of aluminum and its alloys is produced using continuous casting out of a single melt.

Homogeneity

Homogeneity testing is performed by means of spark emission spectroscopy. Tests involve making multiple measurements on individual samples taken at regular intervals along the entire length of each cast rod. Depending on the mass content of the element, the relative standard deviation of multiple measurements between discs or within one disc is typically found between 0.3% - 1% for alloying and other elements and 0.5% - 5% for trace elements.

Analysis

The values listed in this analysis certificate are the results of multiple analyses performed in our chemical analysis laboratory. The analyses are based on established wet chemical procedures.

Description of Sample


This reference material is available in the form of discs (approx. Ø 60 x 25 mm).

Intended use and Stability

This certified reference material is primarily intended for use in spark optical emission spectroscopy. Other applications are X-ray fluorescence spectrometry (XRF) and classical wet chemical procedures. The minimum sample size for wet chemical analysis is 0.2g. The material will remain stable for the period given below (certification validity) if it is stored in a dry and clean environment at room temperature.

Instructions for Use

Calibration measurements should be made within a ring between 2mm and 22mm from the edge of the CRM face. For wet chemical analysis chips have to be prepared by turning or milling of the sample surface.

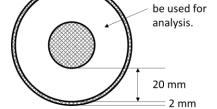

Patrik Bachmann

Head of Inorganic Analytics

White area to

Traceability

Traceability of the certified mass contents to the SI (Système International d'Unités) is ensured by calibration using certified standard solutions or pure metals or substances of known stoichiometry.

Dr. Benedikt Moser CTO

Suisse Technology Partners Ltd. Querstrasse 5 8212 Neuhausen am Rheinfall Switzerland

Phone: +41 52 551 11 00 Fax: +41 52 551 11 99 Email: refmat@suisse-tp.ch

Internet: https://reference-materials.ch

22-Sep-1983 31-Mar-2020 Sep-2058